TECHSPEC[®] ReflX[™] OBJECTIVES #89-725 • 30.56mm WD • 10X/0.23NA UV ENHANCED ALUMINUM

Reflective Objectives are commonly used in applications requiring superior image quality or focusing performance over a broad spectral range. The schwarzschild-style objective exhibits near-diffraction limited performance over the full reflecting range of the chosen coating. By eliminating the problems of chromatic aberration and material absorption associated with standard microscope objectives, these components are ideal for applications requiring high throughput and excellent resolution in the UV or IR, including FTIR spectroscopy, ellipsometry, photolithography, and semiconductor inspection.

Style:	Finite Conjugate
Magnification:	10X
Numerical Aperture NA:	0.23
Working Distance (mm):	30.56
Focal Length FL (mm):	18.4
Field of View, ¾" Sensor:	0.88 x 0.66mm
Field of View, ½" Sensor:	0.64 x 0.48mm
Transmitted Wavefront, RMS:	λ/10
Obscuration (%):	27
Entrance Pupil (mm):	8.4641 (Diameter)

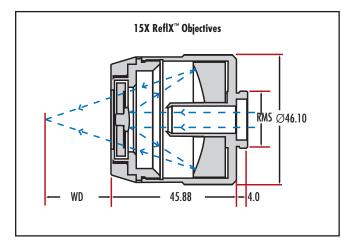


Figure 1: Dimensions of Objective.

Aperture Diameter (mm):	8.6
Diameter of Small Mirror (mm):	8.9
Reflection (%):	0.85
Coating:	UV Enhanced Aluminum
Coating Specification:	R _{ovg} >85% @ 250 - 700nm
Mounting Threads:	RMS
Туре:	Microscope Objective
Wavelength Range (µm):	0.2 - 11
Wavelength Range (nm):	200 - 11000
Manufacturer:	EO

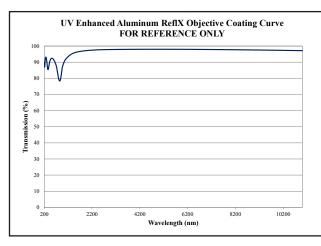


Figure 2: Spectral Reflectance Curve.

